

## DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

### Letermovir

Column: CHIRALPAK® IC-3 [4.6 x 250 mm, 3 µm]

Part No: 83525

Mobile Phase: n-Hexane : Ethanol : Formic acid : DEA (70/30/0.1/0.1,v/v/v/v)


Column Oven Temperature: 25°C

UV Detection: 220nm

Flow Rate: 1.0 ml/min

Diluent: n-Hexane: Ethanol (50/50,v/v)

Sample Concentration: 0.25 mg/ml

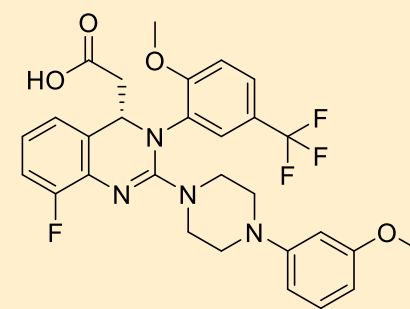

tR<sub>1</sub>: 5.7

N<sub>1</sub>: 6091

T<sub>1</sub>:1.2

tR<sub>2</sub>: 7.2

N<sub>2</sub>: 4083

T<sub>2</sub>:1.2

Rs: 4.1

