

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Isavuconazole

Column: CHIRALPAK® IC-3 [4.6 x 250 mm, 3 µm]

Part No: 83525

Mobile Phase: n-Hexane : Ethanol : Methanol (70/20/10, v/v/v)


Column Oven Temperature: 25°C

UV Detection: 285nm

Flow Rate: 1.0 ml/min

Diluent: Mobile Phase

Sample Concentration: 0.5 mg/ml

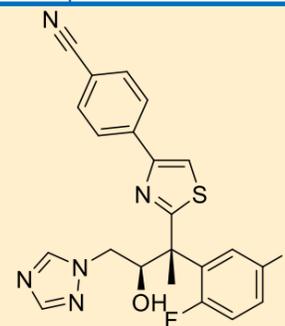

tR₁: 10.9

N₁: 9732

T₁: 1.4

tR₂: 12.6

N₂: 9654

T₂: 1.5

Rs: 3.5

