

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Encorafenib

Column: CHIRALPAK® IH [4.6 x 250 mm, 5 µm]

Part No: 89325

Mobile Phase: n-Hexane/Ethanol/Methanol/Ethanol amine/acetic acid
(80/10/10/0.1/0.1, v/v/v/v/v)


Column Oven Temperature: 40°C

UV Detection: 220nm

Flow Rate: 1.0 ml/min

Diluent: Mobile Phase

Sample Concentration: 0.3 mg/ml

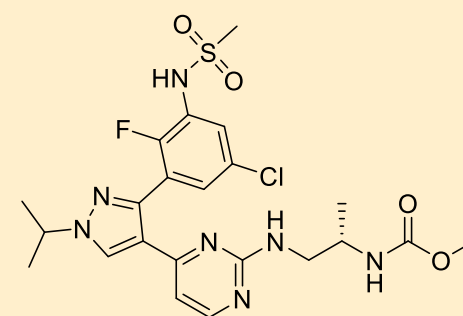

tR ₁ : 9.4	N ₁ : 5389	T ₁ : 1.6	
tR ₂ : 11.5	N ₂ : 4604	T ₂ : 1.5	Rs: 3.5

