

## DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

### Viloxazine

Column: CHIRALPAK® IF [4.6 x 250 mm, 5 µm]

Part No: 86325

Mobile Phase: Methanol: DEA (100/0.1, v/v)


Column Oven Temperature: 25°C

UV Detection: 225nm

Flow Rate: 0.7 ml/min

Diluent: Methanol

Sample Concentration: 1.0 mg/ml

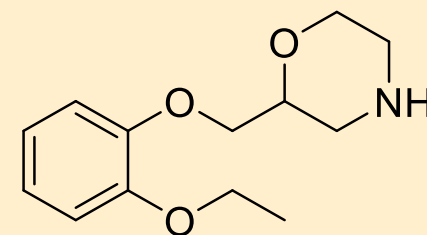

tR<sub>1</sub>: 7.7

N<sub>1</sub>: 8502

T<sub>1</sub>: 1.2

tR<sub>2</sub>: 8.9

N<sub>2</sub>: 7858

T<sub>2</sub>: 1.3

Rs: 3.4

