

## DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

### Dimethindene Maleate

Column: CHIRALPAK® IB N-3 [4.6 x 250 mm, 3 µm]

Part No: 88525

Mobile Phase: Methanol: DEA (100/0.1, v/v)


Column Oven Temperature: 25°C

UV Detection: 260nm

Flow Rate: 1.0 ml/min

Diluent: Mobile Phase

Sample Concentration: 1.0 mg/ml

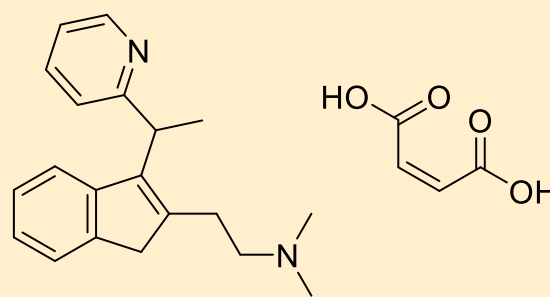

tR <sub>1</sub> : 4.5	N <sub>1</sub> : 8774	T <sub>1</sub> :1.6	
tR <sub>2</sub> : 5.4	N <sub>2</sub> : 10248	T <sub>2</sub> :1.5	Rs: 4.4

