

## DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

### Efinaconazole

Column: CHIRALPAK® IK [4.6 x 250 mm, 5 µm]

Part No: 91325

Mobile Phase: n-Hexane: MTBE: Methanol (90/07/03, v/v/v)


Column Oven Temperature: 15°C

UV Detection: 210nm

Flow Rate: 1.0 ml/min

Diluent: n-Hexane: Ethanol (1:1)

Sample Concentration: 1.0 mg/ml

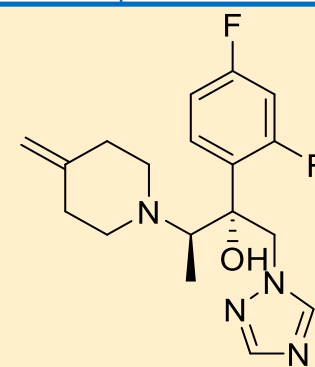

tR <sub>1</sub> : 12.4	N <sub>1</sub> : 11146	T <sub>1</sub> :1.3	
tR <sub>2</sub> : 14.3	N <sub>2</sub> : 10680	T <sub>2</sub> :1.3	Rs: 3.8
tR <sub>3</sub> : 15.7	N <sub>3</sub> : 8638	T <sub>3</sub> :1.1	Rs: 2.1
tR <sub>4</sub> : 17.3	N <sub>4</sub> : 10590	T <sub>4</sub> :1.3	Rs: 2.5

