

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Cabotegravir

Column: CHIRALCEL[®] OX-3R [4.6 x 150 mm, 3 μm]

Part No: 63824

Mobile Phase: Water: ACN: Formic Acid (10/90/0.1,v/v/v)


Column Oven Temperature: 40°C

UV Detection: 260nm

Flow Rate: 0.7 ml/min

Diluent: Mobile Phase

Sample Concentration: 0.5 mg/ml

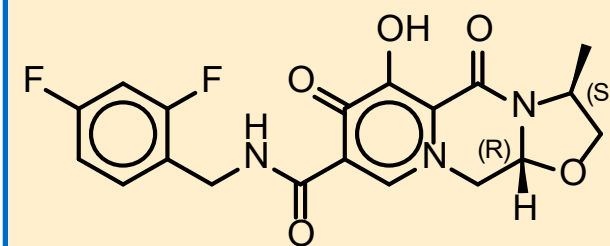

tR ₁ : 7.2	N ₁ : 8069	T ₁ :1.0	
tR ₂ : 9.6	N ₂ : 8546	T ₂ :1.0	Rs: 6.5
tR ₃ : 11.0	N ₃ : 8566	T ₃ :1.0	Rs: 3.1
tR ₄ : 12.3	N ₄ : 8406	T ₄ :1.0	Rs: 2.5

