

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Brivaracetam

Column: CHIRALPAK® IG-U [3.0 x 100 mm, 1.6 µm]

Part No: 87U83

Mobile Phase: 10 mM Ammonium Bicarbonate in Water: ACN (60/40,v/v)


Column Oven Temperature: 25°C

UV Detection: 212 nm

Flow Rate: 0.3 ml/min

Diluent: Mobile Phase

Sample Concentration: 0.5 mg/ml

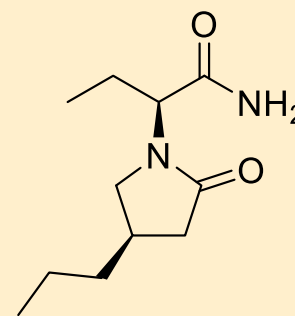

tR ₁ : 3.9	N ₁ : 7455	T ₁ :1.0	
tR ₂ : 4.5	N ₂ :7711	T ₂ :1.0	Rs: 2.6
tR ₃ : 5.8	N ₃ : 8619	T ₃ :1.1	Rs: 5.8
tR ₄ : 8.9	N ₄ : 4375	T ₄ :1.9	Rs: 8.0

