

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Nateglinide

Column: CHIRALPAK® IG-3 [4.6 x 150 mm, 3.0 µm]

Part No: 87524

Mobile Phase A: 0.1% TFA in Water : ACN (50/50, v/v)


Column Oven Temperature: 25°C

UV Detection: 210 nm

Flow Rate: 1.0 ml/min

Diluent: ACN : Water (1:1)

Sample Concentration: 1.0 mg/ml

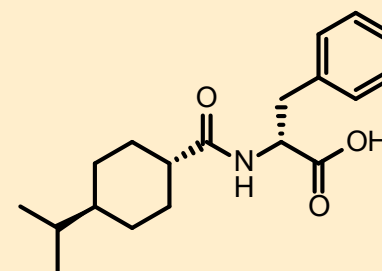

tR ₁ : 5.7	N ₁ : 3634	T ₁ :1.3	
tR ₂ : 8.1	N ₂ :2864	T ₂ :1.2	Rs:4.9

