

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Encorafenib

Column: CHIRALPAK® IH-3 (4.6 x 250mm) 3 μ m Part No: 89525

Mobile Phase: n-Hexane:EtOH:MeOH:Acetic acid:Ethanolamine (80/10/10/0.1/0.1 (v/v/v/v/v))


Column Oven Temperature: 40°C

UV Detection: 220 nm

Flow Rate: 1.0 ml/min

Diluent: n-Hexane: Methanol: Ethanol (80/10/10 (v/v/v))

Sample Concentration: 0.3 mg/ml

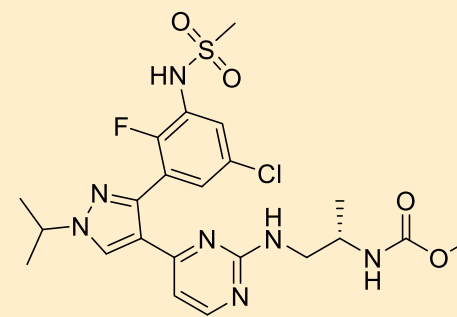

tR₁: 9.9

N₁: 6257

T₁: 1.6

tR₂: 25.6

N₂: 5768

T₂: 1.5

Rs: 3.7

