

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Sotorasib

Column: CHIRALPAK® IK (4.6 X 250mm) 5μ Part No: 91325

Mobile Phase : Acetonitrile


Column Oven Temperature: 25°C

UV Detection: 355 nm

Flow Rate: 1.0 ml/min

Diluent: Acetonitrile

Sample Concentration: 1.0 mg/ml

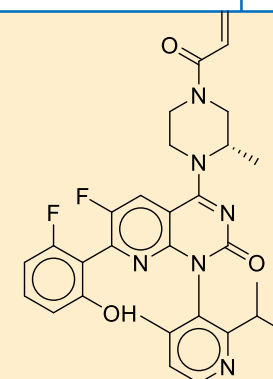

tR_1 : 7.8

N_1 : 4572

T_1 :1.4

tR_2 : 9.8

N_2 :4768

T_2 :1.5

Rs :3.7

