


DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Molnupiravir

Column: CHIRALPAK IC-3 (4.6 X 250mm) 3 μ Part No: 83525
Mobile Phase : n-Hexane/Ethanol/Methanol/Ethanolamine/Acetic acid (60/20/20/0.3/0.3, v/v/v/v/v)
Column Oven Temperature: 40°C
UV Detection: 235 nm
Flow Rate: 1.0 ml/min
Diluent: Mobile Phase
Sample Concentration: 0.5 mg/ml

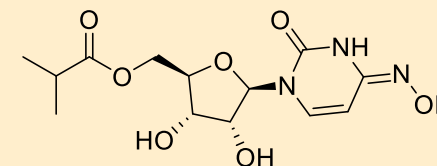

| | | | |
|----------------|--------------|------------|------------|
| t_{R1} : 7.6 | N_1 : 9153 | T_1 :1.3 | |
| t_{R2} : 9.1 | N_2 :10106 | T_2 :1.3 | R_s :4.3 |

