

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Zanubrutinib

Column: CHIRALPAK IB N-5 (4.6 X 250mm) 5 μ Part No: 88325

Mobile Phase : Acetonitrile/Methanol/DEA (80/20/0.1%)


Column Oven Temperature: 25°C

UV Detection: 240 nm

Flow Rate: 1.0 ml/min

Diluent: Mobile phase

Sample Concentration: 1.0 mg/ml

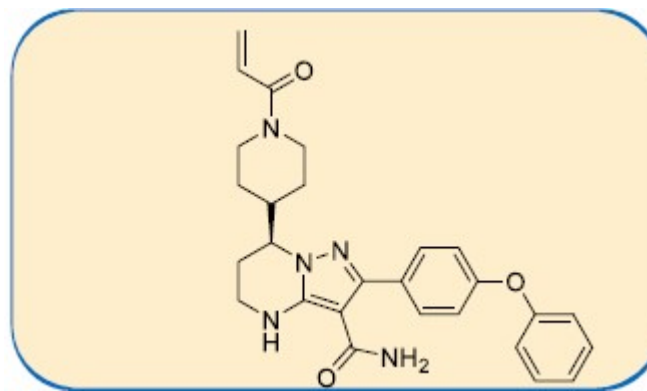

tR₁: 10.5

N₁: 7382

T₁:1.2

tR₂: 11.6

N₂:7292

T₂:1.1

Rs:1.9

