


## DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

### Clomiphene

#### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG [4.6 x 250mm, 5µm] Part No: 87325  
Mobile Phase: Methanol/DEA (100/0.1, v/v)  
Column Oven Temperature: 25°C  
UV Detection: 220 nm  
Flow Rate: 1.0 ml/min  
Diluent: Mobile Phase  
Sample Concentration: 0.5 mg/ml

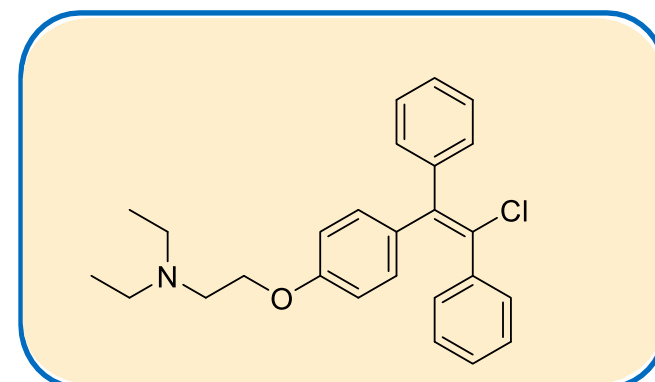

#### SEPARATION PERFORMANCE

$t_{R1}$ : 5.4	$N_1$ :6311	$T_1$ :1.2	
$t_{R2}$ : 6.6	$N_2$ :5494	$T_2$ :1.2	$R_s$ :3.6

