

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Rivaroxaban USP monograph method

Column: CHIRALPAK® IC [4.6 x 250 mm, 5 µm] Part No: 83325

Mobile Phase: Acetonitrile


Column Oven Temperature: 30°C

UV Detection: 254 nm

Flow Rate: 0.7 ml/min

Diluent: Mobile Phase

Sample Concentration: 0.5 mg/ml

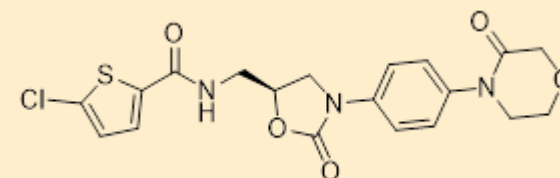

t_{R1} : 21.2

N_1 : 9488

T_1 : 1.1

t_{R2} : 25.8

N_2 : 7122

T_2 : 1.6

R_s : 4.3

