


## DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

### Edoxaban

#### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IK-3 [4.6 x 250mm, 3µm] Part No: 91525  
 Mobile Phase: ACN/Ethanol/2-Butanol/DEA (75/10/15/0.1, v/v/v/v)  
 Column Oven Temperature: 25°C  
 Detection: 290 nm  
 Flow Rate: 0.5 ml/min  
 Diluent: Mobile Phase  
 Sample Concentration: 1.0 mg/ml

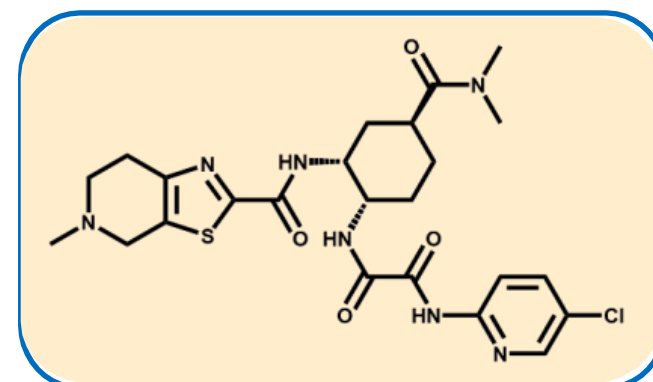

#### SEPARATION PERFORMANCE

t <sub>R1</sub> :10.2	N <sub>1</sub> :15629	T <sub>1</sub> :1.3	
t <sub>R1</sub> :12.2	N <sub>1</sub> :14114	T <sub>1</sub> :1.2	Rs:5.6
t <sub>R1</sub> :12.3	N <sub>1</sub> :15012	T <sub>1</sub> :-	Rs:1.7
t <sub>R1</sub> :13.4	N <sub>1</sub> :15603	T <sub>1</sub> :-	Rs:1.1
t <sub>R1</sub> :14.2	N <sub>1</sub> :16877	T <sub>1</sub> :1.1	Rs:1.7
t <sub>R1</sub> :14.9	N <sub>1</sub> :16534	T <sub>1</sub> :1.1	Rs:1.5
t <sub>R1</sub> :15.8	N <sub>1</sub> :14162	T <sub>1</sub> :1.2	Rs:1.9
t <sub>R1</sub> :24.6	N <sub>1</sub> :17628	T <sub>1</sub> :1.2	Rs:13.7

