


DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Maralixibat

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IN [4.6 x 250mm, 5µm]
Part No: 93325
Mobile Phase: MtBE/MeOH/TFA/Ethanolamine (85/15/0.3/0.3, v/v/v/v)
Column Oven Temperature: 40°C
UV Detection: 220 nm
Flow Rate: 1.0 ml/min
Diluent: Methanol: EtOH (1:4)
Sample Concentration: 1.0 mg/ml

Maralixibat

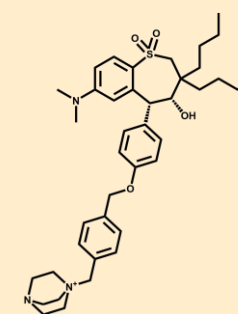

SEPARATION PERFORMANCE

t _{R1} : 5.7	N ₁ : 3385	T ₁ :1.7	
t _{R1} : 7.4	N ₁ : 3250	T ₁ :1.3	Rs: 3.6

