


DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Lisdexamfetamine dimesylate

CHROMATOGRAPHIC CONDITIONS

Column: CROWNPAK® CR-I (+) [3.0 x 150mm, 5 µm]
Part No: 53784
Mobile Phase: Aq. Perchloric acid (pH=1.0): ACN (60/40, v/v)
Column Oven Temperature: 25°C
UV Detection: 215 nm
Flow Rate: 0.5 ml/min
Diluent: Mobile phase
Sample Concentration: 1.0 mg/ml

Lisdexamfetamine dimesylate

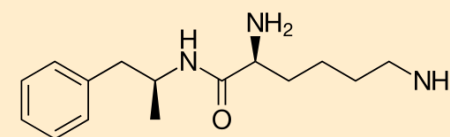

SEPARATION PERFORMANCE

t_{R1} : 2.6	N_1 : 2327	T_1 :1.0	
t_{R1} : 4.3	N_1 : 2536	T_1 :1.1	R_s : 4.8

