


DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Acenocoumarol

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 250mm, 3 µm]
Part No: 87525
Mobile Phase: Methanol: ACN: water: formic acid (50/50/0.2/0.1, v/v/v/v)
Column Oven Temperature: 40°C
UV Detection: 280 nm
Flow Rate: 1.0 ml/min
Diluent: Methanol: ACN (1:1)
Sample Concentration: 1.0 mg/ml

Acenocoumarol

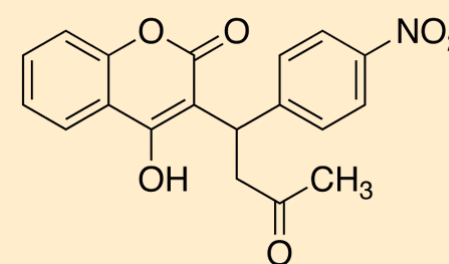

SEPARATION PERFORMANCE

t_{R1} : 4.0	N_1 : 1808	T_1 :1.3	
t_{R1} : 5.1	N_1 : 1615	T_1 :1.2	R_s : 2.6

